

# ELISA Kit for Serum Amyloid A (SAA)

# **Catalog No: tcue298**

**Available Sizes** 

#### **Size:** 96T



**Specifications** 

#### **Research Area:**

Infection immunity;Kidney biomarker;Hepatology;

Species Reactivity:

### Mus musculus (Mouse)

#### Sample Type:

serum, plasma, tissue homogenates, cell lysates, cell culture supernates and other biological fluids

#### Sensitivity:

The minimum detectable dose of this kit is typically less than 0.67ng/mL

#### **Detection Range:**

1.56-100ng/mL

#### **Assay Time:**

3h

#### **Detection Method:**

Enzyme-linked immunosorbent assay for Antigen Detection.

## **Tested Application:**

ELISA

#### SwissProt:


P05366

## **Test Principle**

The test principle applied in this kit is Sandwich enzyme immunoassay. The microtiter plate provided in this kit has been pre-coated



with an antibody specific to Serum Amyloid A (SAA). Standards or samples are then added to the appropriate microtiter plate wells with a biotin-conjugated antibody specific to Serum Amyloid A (SAA). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added, only those wells that contain Serum Amyloid A (SAA), biotin-conjugated antibody and enzyme-conjugated Avidin will exhibit a change in color. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of 450nm  $\pm 10$ nm. The concentration of Serum Amyloid A (SAA) in the samples is then determined by comparing the O.D. of the samples to the standard curve.

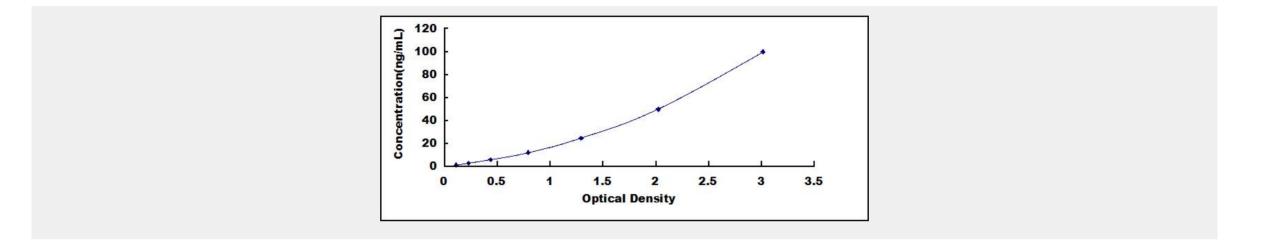

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!

Copyright 2021 Taiclone Biotech Corp.