

# ELISA Kit for Thrombin Activatable Fibrinolysis Inhibitor (TAFI)

# Catalog No: tcue381

Available Sizes

**Size:** 96T



**Research Area:** Metabolic pathway;Hematology;

# **Species Reactivity:**

Homo sapiens (Human)

#### Sample Type:

plasma

# Sensitivity:

The minimum detectable dose of this kit is typically less than 0.133ng/mL

#### **Detection Range:**

0.312-20ng/mL

# **Assay Time:** 3h

### **Detection Method:**

Enzyme-linked immunosorbent assay for Antigen Detection.

# **Tested Application:**

ELISA

# SwissProt:


Q96IY4

Copyright 2021 Taiclone Biotech Corp.



# **Test Principle**

The test principle applied in this kit is Sandwich enzyme immunoassay. The microtiter plate provided in this kit has been pre-coated with an antibody specific to Thrombin Activatable Fibrinolysis Inhibitor (TAFI). Standards or samples are then added to the appropriate microtiter plate wells with a biotin-conjugated antibody specific to Thrombin Activatable Fibrinolysis Inhibitor (TAFI). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added, only those wells that contain Thrombin Activatable Fibrinolysis Inhibitor (TAFI), biotin-conjugated antibody and enzyme-conjugated Avidin will exhibit a change in color. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of 450nm  $\pm$  10nm. The concentration of Thrombin Activatable Fibrinolysis Inhibitor (TAFI) in the samples is then determined by comparing the O.D. of the samples to the standard curve.

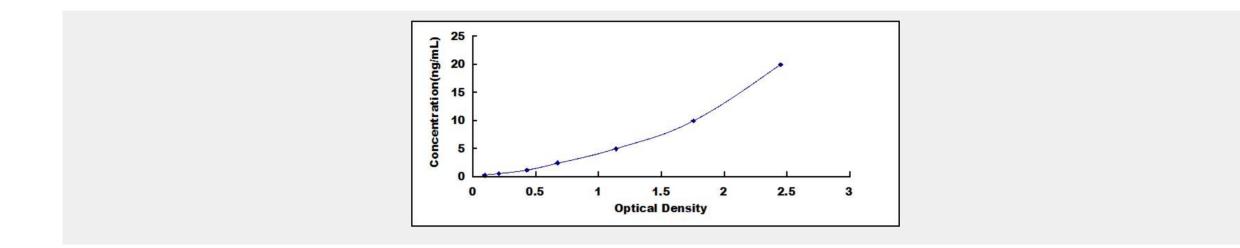

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!

Copyright 2021 Taiclone Biotech Corp.