



# **ELISA Kit for Neural Cell Adhesion Molecule (NCAM)**

**Catalog No: tcue538** 



## **Available Sizes**

Size: 96T



# **Specifications**

#### **Research Area:**

CD & Adhesion molecule; Tumor immunity; Neuro science;

## **Species Reactivity:**

Homo sapiens (Human)

#### **Sample Type:**

Serum, plasma, tissue homogenates and other biological fluids

## **Sensitivity:**

The minimum detectable dose of this kit is typically less than 0.63ng/mL

## **Detection Range:**

1.56-100ng/mL

#### **Assay Time:**

3h

#### **Detection Method:**

Enzyme-linked immunosorbent assay for Antigen Detection.

## **Tested Application:**

**ELISA** 

### **SwissProt:**

P13591


# **Test Principle**

The test principle applied in this kit is Sandwich enzyme immunoassay. The microtiter plate provided in this kit has been pre-coated





with an antibody specific to Neural Cell Adhesion Molecule (NCAM). Standards or samples are then added to the appropriate microtiter plate wells with a biotin-conjugated antibody specific to Neural Cell Adhesion Molecule (NCAM). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added, only those wells that contain Neural Cell Adhesion Molecule (NCAM), biotin-conjugated antibody and enzyme-conjugated Avidin will exhibit a change in color. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of  $450 \text{nm} \pm 10 \text{nm}$ . The concentration of Neural Cell Adhesion Molecule (NCAM) in the samples is then determined by comparing the O.D. of the samples to the standard curve.

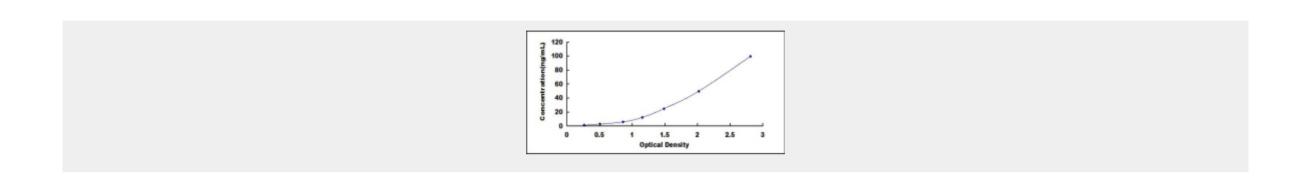

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!