



# **ELISA Kit for Peroxisome Proliferator Activated Receptor Alpha (PPARa)**

Catalog No: tcue771



## **Available Sizes**

Size: 96T



## **Specifications**

#### **Research Area:**

Signal transduction; Tumor immunity; Developmental science;

#### **Species Reactivity:**

Mus musculus (Mouse)

#### **Sample Type:**

Serum, plasma, tissue homogenates, cell lysates, cell culture supernates and other biological fluids

## **Sensitivity:**

The minimum detectable dose of this kit is typically less than 0.056ng/mL

# **Detection Range:**

0.156-10ng/mL

#### **Assay Time:**

3h

#### **Detection Method:**

Enzyme-linked immunosorbent assay for Antigen Detection.

## **Tested Application:**

ELISA

#### **SwissProt:**


P23204





## **Test Principle**

The test principle applied in this kit is Sandwich enzyme immunoassay. The microtiter plate provided in this kit has been pre-coated with an antibody specific to Peroxisome Proliferator Activated Receptor Alpha (PPARa). Standards or samples are then added to the appropriate microtiter plate wells with a biotin-conjugated antibody specific to Peroxisome Proliferator Activated Receptor Alpha (PPARa). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added, only those wells that contain Peroxisome Proliferator Activated Receptor Alpha (PPARa), biotin-conjugated antibody and enzyme-conjugated Avidin will exhibit a change in color. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of  $450 \text{nm} \pm 10 \text{nm}$ . The concentration of Peroxisome Proliferator Activated Receptor Alpha (PPARa) in the samples is then determined by comparing the O.D. of the samples to the standard curve.

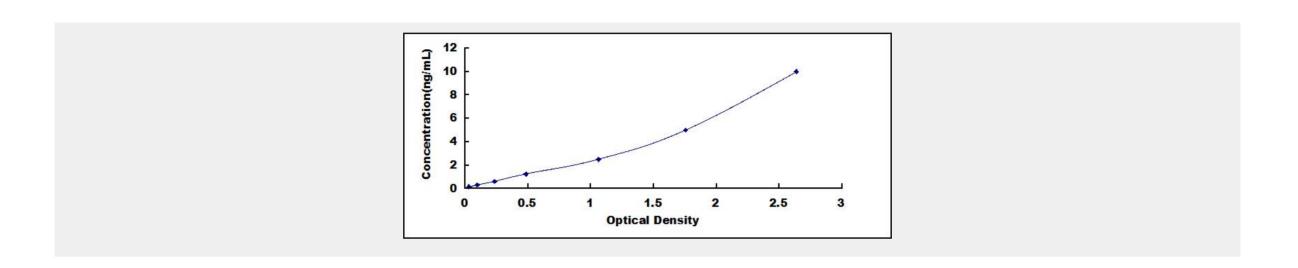

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!