

# ELISA Kit for Toll Like Receptor 2 (TLR2) Catalog No: tcue1756

**Available Sizes** 

## **Size:** 96T



**Specifications** 

### **Research Area:**

Tumor immunity;Infection immunity;Immune molecule;Pulmonology;

# **Species Reactivity:**

Mus musculus (Mouse)

### Sample Type:

tissue homogenates, cell lysates, cell culture supernates and other biological fluids

## Sensitivity:

The minimum detectable dose of this kit is typically less than 0.138ng/mL

# **Detection Range:**

0.312-20ng/mL

# Assay Time:

3h

### **Detection Method:**

Enzyme-linked immunosorbent assay for Antigen Detection.

# **Tested Application:**

ELISA

#### SwissProt:


Q9QUN7

# **Test Principle**

The test principle applied in this kit is Sandwich enzyme immunoassay. The microtiter plate provided in this kit has been pre-coated



with an antibody specific to Toll Like Receptor 2 (TLR2). Standards or samples are then added to the appropriate microtiter plate wells with a biotin-conjugated antibody specific to Toll Like Receptor 2 (TLR2). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added, only those wells that contain Toll Like Receptor 2 (TLR2), biotin-conjugated antibody and enzyme-conjugated Avidin will exhibit a change in color. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of 450nm  $\pm$  10nm. The concentration of Toll Like Receptor 2 (TLR2) in the samples is then determined by comparing the O.D. of the samples to the standard curve.

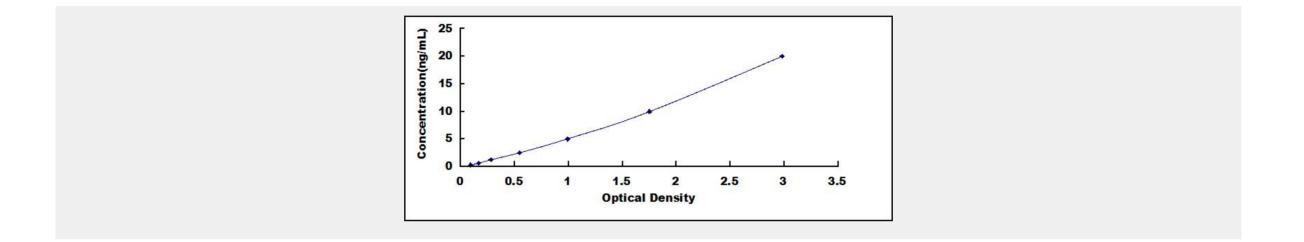

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!

Copyright 2021 Taiclone Biotech Corp.