

ELISA Kit for p53 Upregulated Modulator Of Apoptosis (PUMA)

Catalog No: **tcue1766**



Available Sizes

Size: 96T



Specifications

Research Area:

Apoptosis; Tumor immunity;

Species Reactivity:

Homo sapiens (Human)

Sample Type:

Tissue homogenates, cell lysates and other biological fluids.

Sensitivity:

The minimum detectable dose of this kit is typically less than 0.059ng/mL

Detection Range:

0.156-10ng/mL

Assay Time:

3h

Detection Method:

Enzyme-linked immunosorbent assay for Antigen Detection.

Tested Application:


ELISA

SwissProt:

Q9BXH1

Test Principle

The test principle applied in this kit is Sandwich enzyme immunoassay. The microtiter plate provided in this kit has been pre-coated with an antibody specific to p53 Upregulated Modulator Of Apoptosis (PUMA). Standards or samples are then added to the appropriate microtiter plate wells with a biotin-conjugated antibody specific to p53 Upregulated Modulator Of Apoptosis (PUMA). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added, only those wells that contain p53 Upregulated Modulator Of Apoptosis (PUMA), biotin-conjugated antibody and enzyme-conjugated Avidin will exhibit a change in color. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of $450\text{nm} \pm 10\text{nm}$. The concentration of p53 Upregulated Modulator Of Apoptosis (PUMA) in the samples is then determined by comparing the O.D. of the samples to the standard curve.

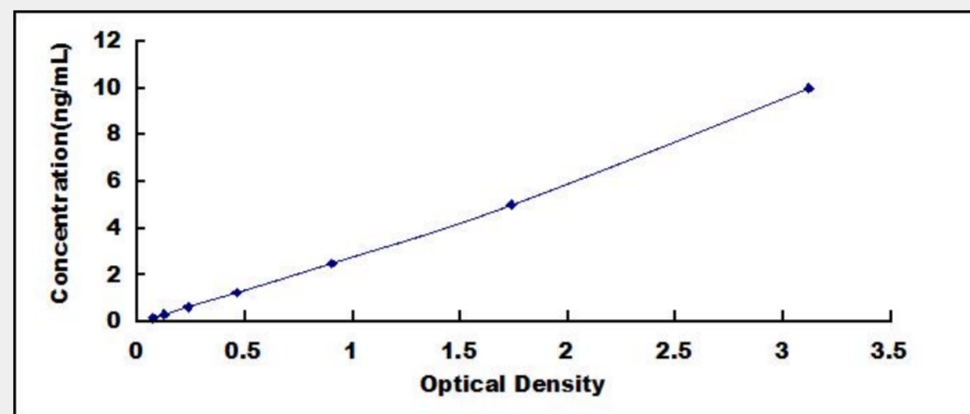

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!