

# ELISA Kit for Peptidase Inhibitor 3, Skin Derived (PI3) Catalog No: tcue2662

**Available Sizes** 

**Size:** 96T

Specifications

**Research Area:** 

Immune molecule;

Species Reactivity: Homo sapiens (Human)

### Sample Type:

Serum, plasma, urine, saliva, tissue homogenates, cell lysates, cell culture supernates and other biological fluids

### Sensitivity:

The minimum detectable dose of this kit is typically less than 0.124ng/mL

### **Detection Range:**

0.312-20ng/mL

### Assay Time:

3h

#### **Detection Method:**

Enzyme-linked immunosorbent assay for Antigen Detection.

## **Tested Application:**

ELISA

#### SwissProt:


P19957

# **Test Principle**

The test principle applied in this kit is Sandwich enzyme immunoassay. The microtiter plate provided in this kit has been pre-coated



with an antibody specific to Peptidase Inhibitor 3, Skin Derived (PI3). Standards or samples are then added to the appropriate microtiter plate wells with a biotin-conjugated antibody specific to Peptidase Inhibitor 3, Skin Derived (PI3). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added, only those wells that contain Peptidase Inhibitor 3, Skin Derived (PI3), biotin-conjugated antibody and enzyme-conjugated Avidin will exhibit a change in color. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of 450nm  $\pm 10$ nm. The concentration of Peptidase Inhibitor 3, Skin Derived (PI3) in the samples is then determined by comparing the O.D. of the samples to the standard curve.

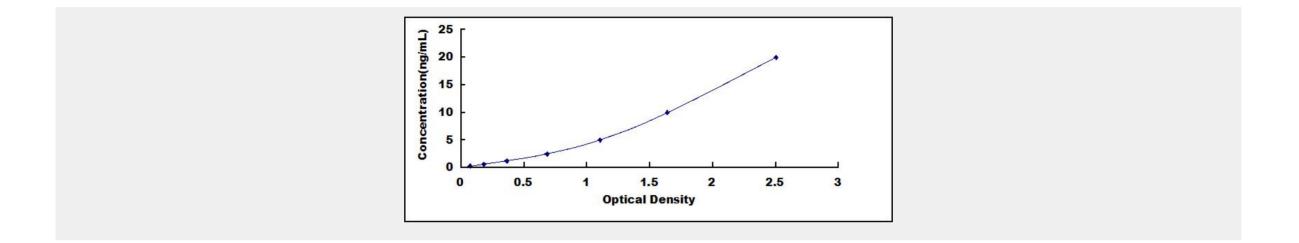

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!

Copyright 2021 Taiclone Biotech Corp.