



ELISA Kit for Wingless Type MMTV Integration Site Family, Member 1 (WNT1)

Catalog No: tcue2939


Av	ailable Sizes
Size: 96T	
Sp	ecifications
Research Tumor imr	
_	Reactivity: iens (Human)
Sample T	Type: nogenates, cell lysates and other biological fluids
Sensitivit The minim	ty: num detectable dose of this kit is typically less than 0.057ng/mL
Detection 0.156-10n	
Assay Tir 3h	ne:
	n Method: nked immunosorbent assay for Antigen Detection.
Tested A ELISA	pplication:
SwissPro	t:





Test Principle

The test principle applied in this kit is Sandwich enzyme immunoassay. The microtiter plate provided in this kit has been pre-coated with an antibody specific to Wingless Type MMTV Integration Site Family, Member 1 (WNT1). Standards or samples are then added to the appropriate microtiter plate wells with a biotin-conjugated antibody specific to Wingless Type MMTV Integration Site Family, Member 1 (WNT1). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added, only those wells that contain Wingless Type MMTV Integration Site Family, Member 1 (WNT1), biotin-conjugated antibody and enzyme-conjugated Avidin will exhibit a change in color. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of 450nm \pm 10nm. The concentration of Wingless Type MMTV Integration Site Family, Member 1 (WNT1) in the samples is then determined by comparing the O.D. of the samples to the standard curve.

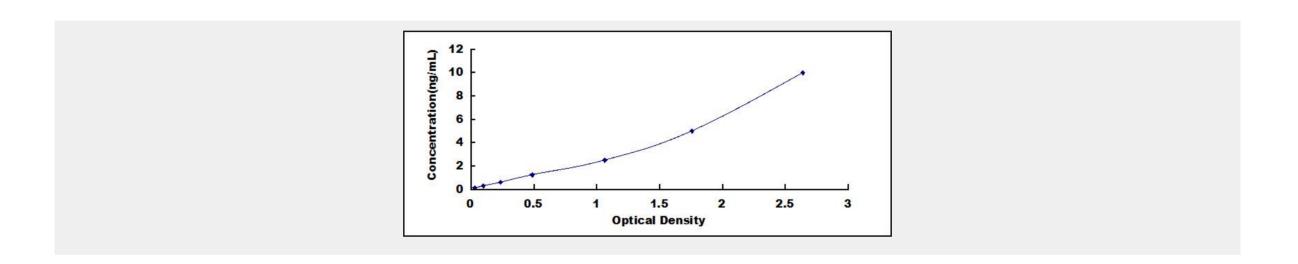

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!