



ELISA Kit for A Disintegrin And Metalloprotease 22 (ADAM22)

Catalog No: tcue3435


| Available Sizes |
|---|
| ize: 96T |
| Specifications |
| esearch Area: nzyme & Kinase; |
| pecies Reactivity: omo sapiens (Human) |
| ample Type: erum, plasma, tissue homogenates and other biological fluids |
| ensitivity: he minimum detectable dose of this kit is typically less than 0.058ng/mL |
| Petection Range: .156-10ng/mL |
| ssay Time: h |
| Petection Method: nzyme-linked immunosorbent assay for Antigen Detection. |
| ested Application: LISA |
| wissProt: |





Test Principle

The test principle applied in this kit is Sandwich enzyme immunoassay. The microtiter plate provided in this kit has been pre-coated with an antibody specific to A Disintegrin And Metalloprotease 22 (ADAM22). Standards or samples are then added to the appropriate microtiter plate wells with a biotin-conjugated antibody specific to A Disintegrin And Metalloprotease 22 (ADAM22). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added, only those wells that contain A Disintegrin And Metalloprotease 22 (ADAM22), biotin-conjugated antibody and enzyme-conjugated Avidin will exhibit a change in color. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of 450nm \pm 10nm. The concentration of A Disintegrin And Metalloprotease 22 (ADAM22) in the samples is then determined by comparing the O.D. of the samples to the standard curve.

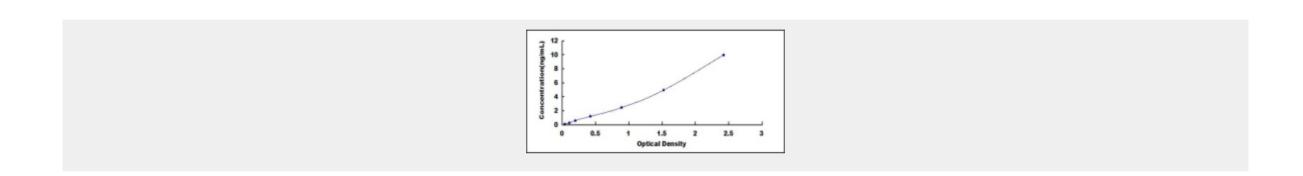

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!