



ELISA Kit for Metallothionein 1E (MT1E)

Catalog No: tcue3444

Avai	ilable Sizes
Size: 96T	
Spec	cifications
Research A Tumor immu	
Species Rea Homo sapier	
Sample Typ Serum, plasr	pe: ma, tissue homogenates and other biological fluids.
Sensitivity: The minimur	: m detectable dose of this kit is typically less than 0.112ng/mL
Detection F 0.312-20ng/i	
Assay Time 3h	2:
Detection N Enzyme-linke	Method: ed immunosorbent assay for Antigen Detection.
Tested App ELISA	olication:
SwissProt:	

Test Principle


P04732

The test principle applied in this kit is Sandwich enzyme immunoassay. The microtiter plate provided in this kit has been pre-coated





with an antibody specific to Metallothionein 1E (MT1E). Standards or samples are then added to the appropriate microtiter plate wells with a biotin-conjugated antibody specific to Metallothionein 1E (MT1E). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added, only those wells that contain Metallothionein 1E (MT1E), biotin-conjugated antibody and enzyme-conjugated Avidin will exhibit a change in color. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of $450 \, \text{nm} \pm 10 \, \text{nm}$. The concentration of Metallothionein 1E (MT1E) in the samples is then determined by comparing the O.D. of the samples to the standard curve.

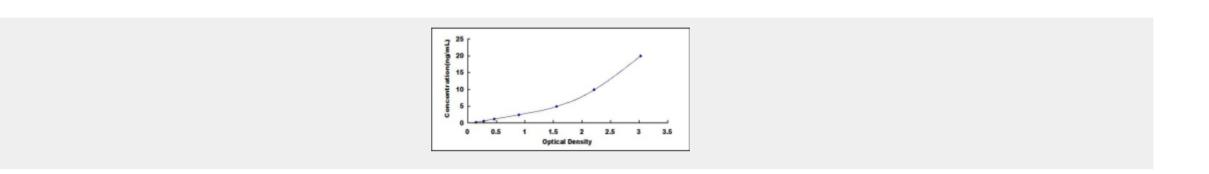

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!