

ClonFluor[®] Mounting Medium

Catalog No: tcst2080



Available Sizes

Size: 25ml



Specifications

Storage Buffer:

Contain NaN₃ as a preservative

Tested Application:

IHC-FS;IHC-Whole Mounting;(F)ISH;ICC/IF

Storage Instruction:


Avoid exposure to light, store at Room Temperature. (22-25°C)

Protocol:

Prepare samples and perform staining according to established procedures Following aqueous final step, remove excess liquid around the sample and add 1 drop (~50 µL) of DAPI ClonFluor[®] Mounting Media for an area of 22 x 22 mm Mount coverslip and remove.

Product Description

ClonFluor[®] Mounting Medium is a water-soluble compound recommended for slides mounted after a staining procedure having an aqueous final step. Since ClonFluor[®] Mounting Medium is water-soluble, the coverslip may be removed by submerging the slide in a PBS solution until the coverslip is loosened. This mounting medium also provides a semi-permanent seal for storage of slide preparations.

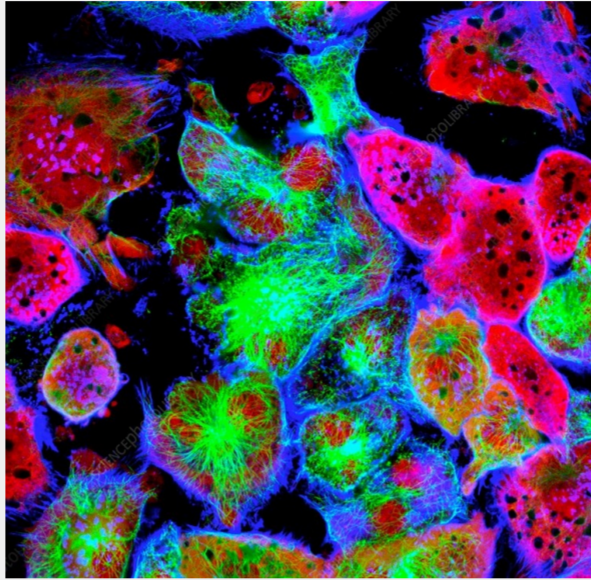

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!