



**CHO-K1** 

**Catalog No: tcel62** 



**Available Sizes** 

Size: 1×10<sup>6</sup>cells/t25culturebottle



## **Specifications**

## **Application:**


transfection host

## **Subculturing:**

Remove and discard culture medium. Briefly rinse the cell layer with DPBS solution to remove all traces of serum that contains trypsin inhibitor. Add 1.0 to 2.0 mL of Trypsin-EDTA solution to flask and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes). Cells that are difficult to detach may be placed at 37°C to facilitate dispersal. Add 4.0 to 6.0 mL of complete growth medium and aspirate cells by gently pipetting. Add appropriate aliquots of the cell suspension to new culture vessels.

## **Product Description**


The cells require Proline in the medium for growth. The CHO-K1 cell line was derived as a subclone from the parental CHO cell line initiated from a biopsy of an ovary of an adult Chinese hamster by T. T. Puck in 1957. The cells require Proline in the medium for growth.

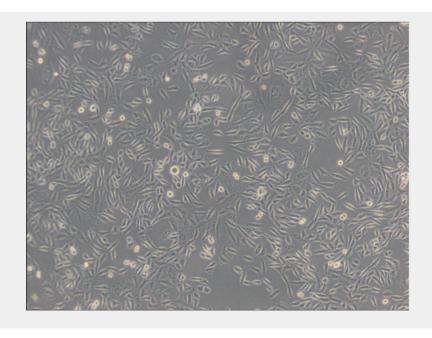










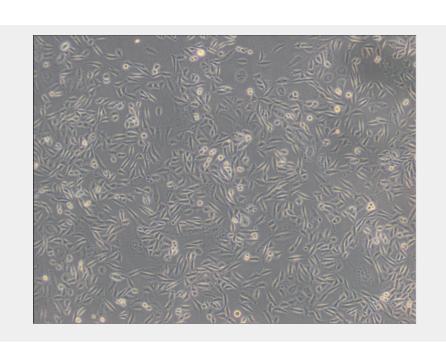

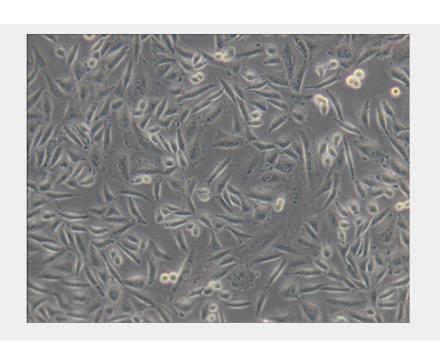

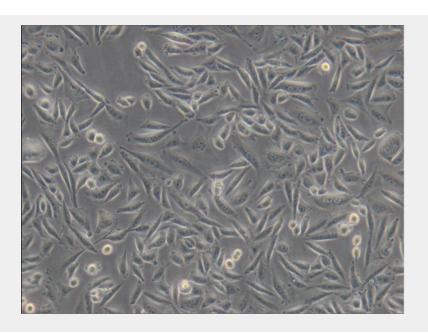

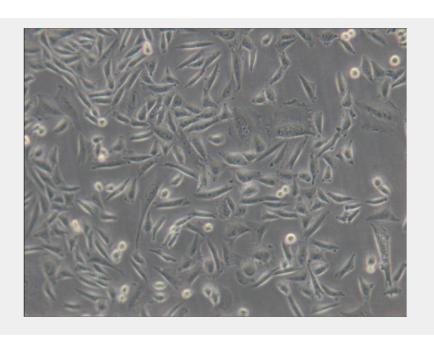

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!