

RL95-2

Catalog No: tcel197



Available Sizes

Size: 1×10⁶ cells/t25 culture bottle




Specifications


Subculturing:

Remove and discard culture medium. Briefly rinse the cell layer with DPBS solution to remove all traces of serum that contains trypsin inhibitor. Add 1.0 to 2.0 mL of Trypsin-EDTA solution to flask and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes). Cells that are difficult to detach may be placed at 37°C to facilitate dispersal. Add 4.0 to 6.0 mL of complete growth medium and aspirate cells by gently pipetting. Add appropriate aliquots of the cell suspension to new culture vessels.

Product Description

The cells possess alpha keratin, well defined junctional complexes, tonofilaments and surface microvilli.

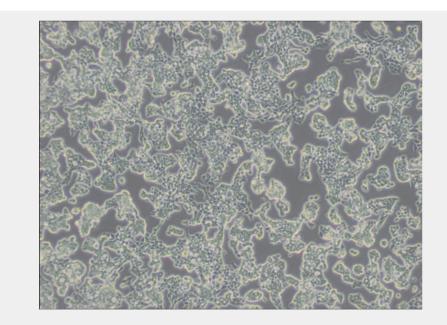


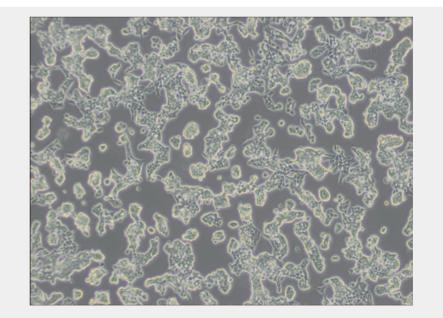

Copyright 2021 Taiclone Biotech Corp.



Web: www.taiclone.com Tel: +886-2-2735-9682 Email: order@taiclone.com

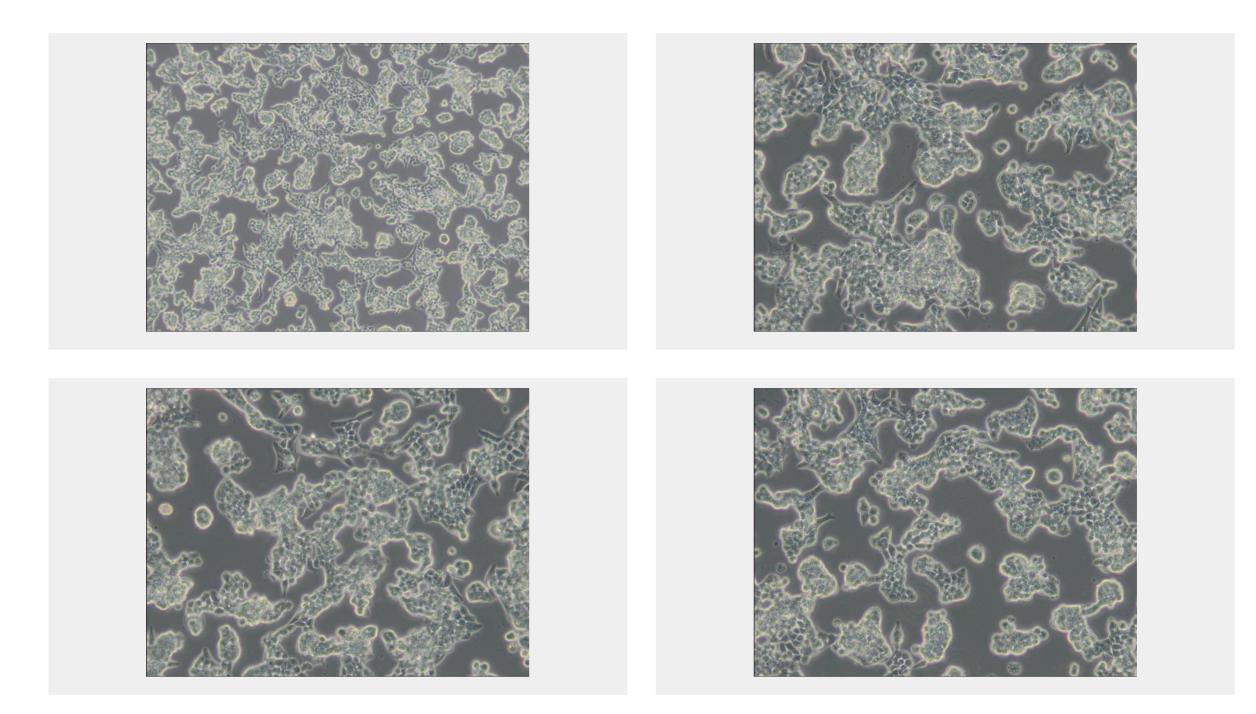

All products are for RESEARCH USE ONLY. Not for diagnostic & therapeutic purposes!

Copyright 2021 Taiclone Biotech Corp.